Warmer Pad

ALMOHADILLA ELÉCTRICA

MANUAL DE INSTRUCCIONES

HEATING PAD

USERS MANUAL

COUSSIN CHAUFFANT

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ALMOFADA ELECTRICA

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CUSCINETTO RISCALDANTE

MANUALE D'ISTRUZIONI

ELEKTRISCHES HEIZKISSEN

ANLEITUNG

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Please read this manual carefully before use,
and keep it for any future reference

TECHNICAL DATA

Model:	Warmer Pad	Power supply:	220-240 V ~ 50Hz
Max. power:	100 W	Dimensions:	40x30 cm

This instruction manual belongs to this device. It contains important information about starting up and operation. Read the instruction manual thoroughly. Non-observance of these instructions can result in serious injury or damage to the device.


Do not puncture the warmer pad!

Do not chlorinate!

Do not use the warmer pad when it is folded!

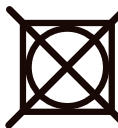Do not dry the warmer pad in a tumble drier! The pad cover may not be put in a dryer!

Not to be used by very young children (0-3 years)!

Do not iron the warmer pad!
You may iron the pad cover at a low temperature!

Only use the warmer pad indoors!

The heating pad can be washed by hand! You may wash the cushion cover in a normal wash at max. 30 °C!

Do not dry clean!

SAFETY PRECAUTIONS

- Carefully check the heat pad before each use, if you notice wear or damage.
- Do not use if you notice wear, damage or signs of improper use on the device, the switch or the cables or if the device is not working.
- Before connecting the device to your power supply, please ensure that the supply voltage stated on the rating plate is compatible with your mains supply.
- Do not fold or crease the heat pad during use.
- Neither pierce the pad with, nor attach safety pins or other spiky or sharp objects to the pad.
- Do not use the device on children or persons who are disabled, asleep or insensitive to heat (persons who cannot react on overheating).
- Children below the age of 3 years are not allowed to use this device, because they are not able to react on overheating.
- The device must not be used by young children over the age of 3 years, unless the control unit has been adjusted accordingly by its parents or legal guardians or unless the child has been well instructed how to use the control unit properly.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- Never use the pad unsupervised.
- The device is only for private use and is not to be used in hospitals.
- Do not fall asleep while the heat pad is switched on.
- Prolonged use of the warmer pad at a high setting could lead to skin burns.
- Do not cover the pad with any other cushions. The control switch should not be covered or placed on or below the pad while the unit is in operation.
- Do not sit on the heat pad. Instead, rest the pad on or against the part of the body to be treated. The intended use of the is to treat only back and neck.
- Never touch a heat pad that has fallen into water. Unplug the unit from the mains outlet immediately.
- Keep the mains lead away from hot surfaces.
- Never carry, pull or turn the heat pad by the mains lead and do not let the lead

get tangled.

- Do not use the heat pad when wet and only use it in dry environments (not in the bathroom or similar environments).
- The controls and leads must not be exposed to any kind of moisture.
- If a fault occurs, do not attempt to repair the unit yourself. Repairs must only be carried out by an authorized specialist dealer or other suitably qualified personnel.
- If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
- A damaged power supply cable must be replaced with a power cable of the same type.
- Heat pad may only be operated with the control unit (SS02) supplied.
- When storing the heat pad, allow it to cool down before folding.
- Should you have any concerns about health issues, consult your doctor before using the heat pad.
- Do not use the heat pad on parts of the body that are swollen, injured or inflamed.
- Should you experience persistent pains in the muscles and joints, please inform your doctor. Persistent pains can be a symptom of serious illness.
- If you experience any pain or discomfort while using the unit, stop using it immediately.

OPERATION

BEFORE THE FIRST USE

Inspect the product completely for damage. If the contents are incomplete or damaged, contact your retailer immediately.

USE

1. Connect the power plug to the power outlet.
2. Turn on the product.

HEAT SETTINGS

0 = off

1 = low level

2 = medium level

3 = high level

EN

- It takes approximately 15 minutes for the pad to reach the desired temperature. Remember that a temperature that feels hot in your hand can be extremely hot in other parts of your body.
- You may notice a slight odor the first time the pad is heated. This smell is completely harmless, but be sure to air the room after the blanket has been heated for the first time.
- Start by heating the pad to the highest heat level. After 15 minutes the temperature can be adjusted to the desired level.
- The heating function of the pad turns off automatically after approximately 90min. If you wish to reheat the pad after approximately 90min, simply set the heat level back to O for approximately 5 seconds. You can then select the desired heat level again.

CLEANING & MAINTENANCE

- Before cleaning the pad, unplug it from the mains outlet and let it cool down for at least ten minutes.
- The heating pad is equipped with a detachable cable. Disconnect the plug and remove the cable from the heating pad.
- You may wash the heating pad gently by hand. The best way to do this is to put the heating pad in a bath with lukewarm water and some mild detergent and then gently squeeze it out.
- Rinse the heating pad several times to remove all the detergent.
- Store the unit in a clean and dry place, laid out flat with nothing covering or resting on top of it.
- Only use the appliance once it is completely dry.

WARRANTY

The benefits conferred by this warranty are in addition to all implied warranties, other rights and remedies in respect of the product which the consumer has under the Trade Practices Act and similar State and Territory Laws. The original purchaser of this product is provided with the follow warranty, subject to the following conditions. This product is warranted for a period of 2 years from faulty, the product will need to be returned to the retailer where it was purchased to be replaced free of charge or a refund issued. In order to obtain an exchange or a refund, proof of purchase must be presented to the retailer.

DISCARDING WORN HOUSEHOLD DEVICES

The European directive 2012/19/EC pertaining to electrical and electronic waste prohibits such household goods from being discarded with local refuse services. Worn devices must be collected separately to optimize the recovery and recycling rates of the components therein to reduce detrimental effects on the environment and on human health.

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